



Email : sales@scilabexport.com

Phone: +91-7082934803


Product Name :

Optical Lens Set

Product Code :

SCL255 /042

Image not found or type unknown

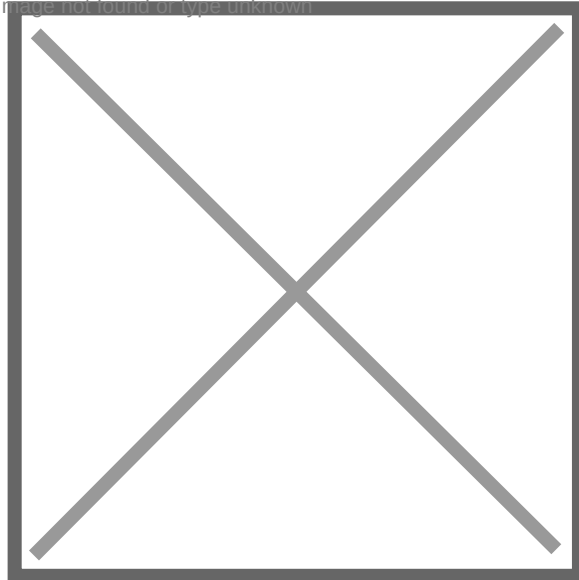

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market | 133001, Harvana



Email : sales@scilabexport.com

Phone: +91-7082934803

Description :

Optical Lens Set (Product Code: SCL255 /042) is a 225-piece trial lens set used to check refraction and try corrective lens powers in eye clinics, hospital eye departments and optometry colleges. It combines 32 pairs each of convex and concave spherical lenses, 20 pairs each of convex and concave cylindrical lenses, a prism series and a trial frame, in glass and metal. It is supplied by **Jain Scientific Equipments Pvt Ltd (ScienceLab)**, Ambala, India.

Specifications

Parameter	Detail
Total pieces	225
Material	Glass and metal
Rim option	Adjustable optometric metal rim cover available
Spherical lenses	32 pairs each, convex (cx) and concave (cc)
Spherical powers	0.25 to 4.00 in 0.25 steps; 4.50 to 6.00 in 0.50 steps; 7.00 to 14.00 in 1.00 steps; 16.00 to 20.00 in 2.00 steps
Cylindrical lenses	20 pairs each, convex (cx) and concave (cc)
Cylindrical powers	0.25 to 3.50 in 0.25 steps; 4.00 to 6.00 in 0.50 steps
Prisms	Listed values 1/2, 1, 2, 3, 4, 5, 6, 8, 10 and 12
Accessories	Maddox rod, red glass, green glass, stenopaic slit, pinhole, 2 blank discs
Frame	Trial frame included
Documentation	User manual in English
Product code	SCL255 /042
Supplied by	Jain Scientific Equipments Pvt Ltd (ScienceLab), Ambala, India

What Is a Trial Lens Set Used For?

In subjective refraction the examiner places spherical lenses in the trial frame and changes them step by step until the patient reports the clearest vision, then adds cylindrical lenses to refine the correction for astigmatism. Lens powers in trial sets are read in dioptres, and the 0.25 steps at the low end of this set allow fine adjustments. The prisms are used when checking how the two eyes work together.

The accessories widen the range of tests: the Maddox rod and the red and green glasses for binocular balance checks, the stenopaic slit for locating the axis of astigmatism, and the pinhole for judging whether reduced vision is optical in origin. Optometry colleges use the Optical Lens Set to train students in these routines. Hold lenses by the rim and return each one to its own slot so the set stays complete.

What to

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

- The

1226/1-5, Bengali Mohalla, Science Market | 133001, Haryana



Email : sales@scilabexport.com

Phone: +91-7082934803