


**Product Name :**

Electric Field Apparatus

**Product Code :**

SLE/EAM/15

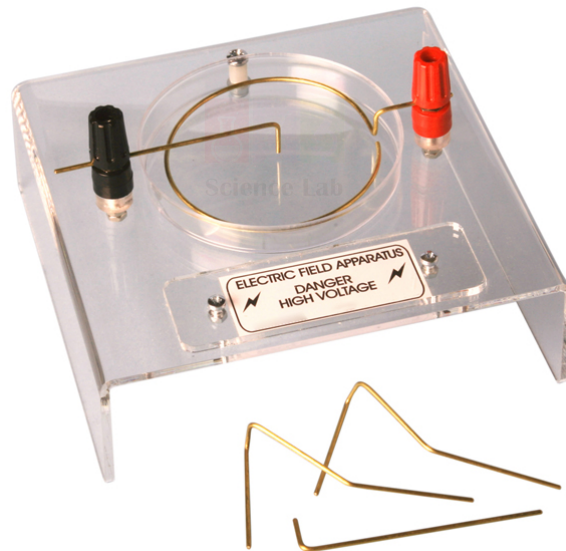

**Description :**

Place liquid paraffin in the 90mm petri dish and sprinkle semolina or fine grass seed (not supplied) on the surface.

A simple model for demonstrating electric fields between various electrode configurations.

Apply a high voltage via the 4mm terminals and the field pattern is generated.

The clear acrylic base allows the pattern to be projected via an overhead projector. (150 x 150 x 60).

Note: This apparatus requires a safe, low-current limited 0-5kV EHT supply. Investigate electric fields

Contact Science Lab for your Educational School Science Lab Equipments. We are best educational lab equipments exporter, educational lab equipments manufacturers, educational lab equipments supplier, educational lab equipments suppliers, educational lab export, educational lab exporter.