



Email : sales@scilabexport.com

Phone: +91-7082934803


Product Name :

Gaspak EZ Anaerobe Sachet with Indicator

Product Code :

SCI-241/039

Image not found or type unknown

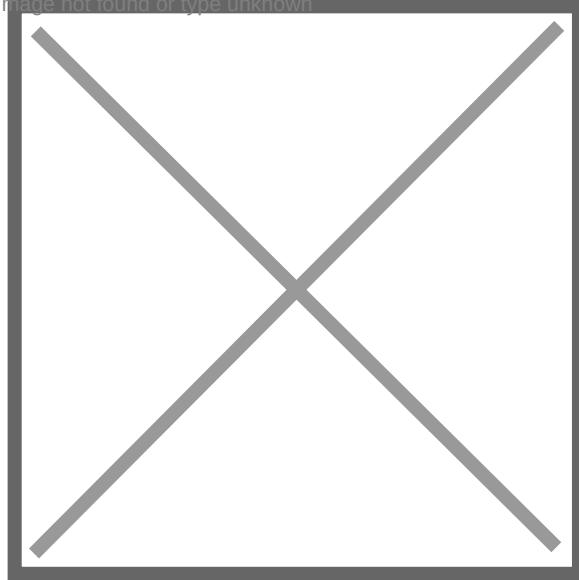

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market | 133001, Haryana



Email : sales@scilabexport.com

Phone: +91-7082934803

Description :

Gaspak EZ Anaerobe Sachet with Indicator (Product Code: SCI-241/039) is a gas-generating sachet used to create an anaerobic atmosphere inside the Gas Pak 100 anaerobic system when culturing anaerobic bacteria. Placed in the sealed jar, it produces anaerobic conditions within about 2½ hours and more than 15% carbon dioxide within 24 hours, and an oxygen indicator is included so staff can check the result. It is supplied by **Jain Scientific Equipments Pvt Ltd (ScienceLab)**, Ambala, India.

Specifications

Parameter	Detail
For use with	Gas Pak 100 anaerobic system
Function	Creates an anaerobic atmosphere in the sealed jar
Time to anaerobic atmosphere	Within about 2½ hours
Carbon dioxide	Greater than 15% within 24 hours
Indicator	Oxygen indicator included
Shelf life	Shelf life depends on the product and its storage conditions; ask us for the expiry date of the batch you are quoted.
Product code	SCI-241/039
Supplied by	Jain Scientific Equipments Pvt Ltd (ScienceLab), Ambala, India

What Is the Gaspak EZ Anaerobe Sachet Used For?

The sachet is the active consumable of the anaerobic jar. Sealed in the jar with inoculated plates or tubes, it takes up the oxygen in the enclosed air and releases carbon dioxide, giving obligate anaerobes the conditions they need and providing the carbon dioxide that many fastidious organisms prefer. The oxygen indicator gives a visual check that the atmosphere has become anaerobic before the incubation is relied upon.

Load all plates first, then open the sachet, put it in the jar and close the lid straight away, since the reaction begins once the sachet meets air. A sachet serves one jar run and cannot be regenerated, so college microbiology practicals usually plan sachet stock by counting the anaerobic incubations in the term. Store unopened sachets as directed on the pack.

What to Specify When Ordering

- Number of sachets, based on one per jar run.
- The jar you use, as the sachet is listed for the Gas Pak 100 anaerobic system.
- Whether you need the batch expiry date before ordering stock for a full term.

Supplier

ScienceLab

supplies to

and dealer

or together

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market | 133001, Haryana

aryana,

leges

stocking

on the



Email : sales@scilabexport.com

Phone: +91-7082934803