



Email : sales@scilabexport.com

Phone: +91-7082934803


Product Name :

Funnel

Product Code :

SCL-523/011

Image not found or type unknown

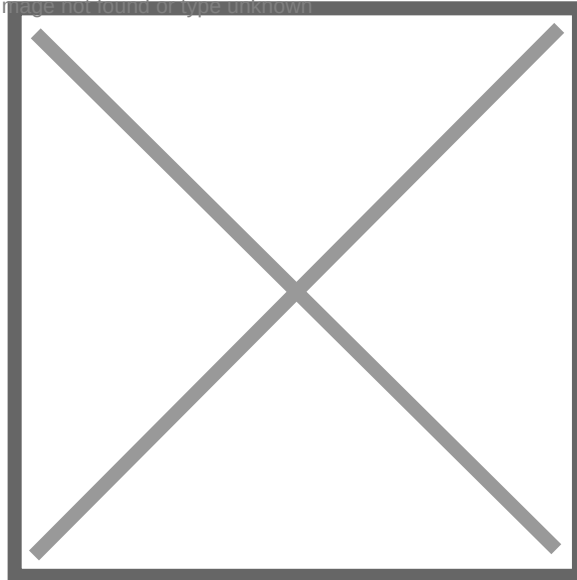

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India



Email : sales@scilabexport.com

Phone: +91-7082934803

Description :

Funnel (Product Code: SCL-523/011) is a long-stem laboratory funnel moulded in clear polypropylene, used for gravity filtration and for pouring liquids into narrow-necked flasks without spills in school and college labs. Its 60° cone, smooth walls and external ribs that prevent air lock keep liquid moving, and it can be steam autoclaved. It is made by **Jain Scientific Equipments Pvt Ltd (ScienceLab)**, an OEM manufacturer, at its own factory in Ambala, India.

Specifications

Parameter	Detail
Material	Plastic — moulded polypropylene, clear
Stem	Long stem
Cone angle	60°
Walls	Smooth
Ribs	External ribs to prevent air lock
Sterilisation	Can be steam autoclaved
Product code	SCL-523/011
Manufacturer	Jain Scientific Equipments Pvt Ltd (ScienceLab) — OEM
Country of origin	India (made in Ambala, Haryana)

What Is the Polypropylene Funnel Used For?

For filtration, a circle of filter paper is folded into a cone, opened out and seated in the funnel; the 60° cone lets the paper lie against the smooth wall so the liquid passes through the paper rather than running down the sides. The long stem fills with filtrate as the work proceeds, which helps draw the liquid through, and the external ribs keep a small air gap in the neck of a conical flask so the receiver can vent. The same funnel is handy for refilling reagent bottles and for pouring solutions into volumetric flasks during titration work.

Being polypropylene, it survives the knocks of a busy practical room far better than glass, which is why schools add it to their lab plastic ware alongside beakers and bottles, but it should not be placed over a flame or on a hot plate. Rinse it straight after use, clean with warm water, a mild detergent and a soft brush, and avoid abrasive powders. Where sterile handling is needed, for example in microbiology classes, steam autoclave it before use.

OEM Manufacturer

Jain Scientific Equipments Private Limited, which trades as ScienceLab, is the OEM (original equipment manufacturer) of this polypropylene Funnel. It is produced at the company's own licensed factory at Saha, Ambala, Haryana, India, and supplied factory-direct to schools, colleges, dealers and importers rather than through a trading intermediary.

Repe
funne

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

D... 1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India



Email : sales@scilabexport.com

Phone: +91-7082934803

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India