



Email : sales@scilabexport.com

Phone: +91-7082934803


Product Name :

Glazed Tiles

Product Code :

SCL-581/040

Image not found or type unknown

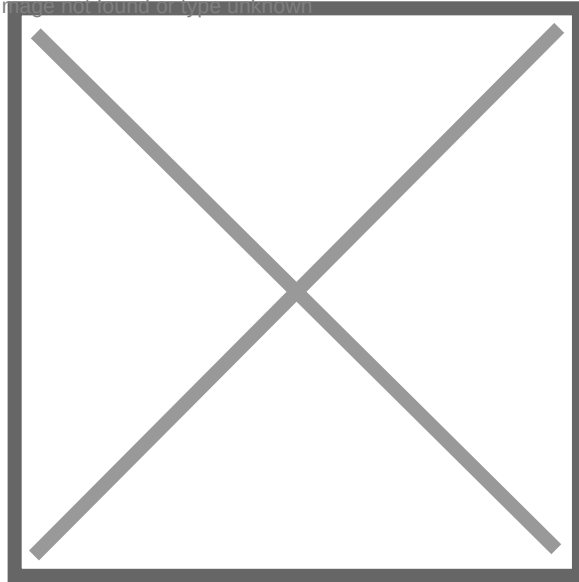

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India



Email : sales@scilabexport.com

Phone: +91-7082934803

Description :

Glazed Tiles (Product Code: SCL-581/040) are white, square ceramic tiles with a crisp, smooth glazed surface that is easy to clean and maintain. In the chemistry laboratory they serve as a clean white background, most often placed under the conical flask during a titration so that the colour change at the end point is easy to judge, and also as a smooth surface for spotting dropsof indicator. A tile for every student station is a small but useful part of a titration kit for school practicals. The tiles are supplied by **Jain Scientific Equipments Pvt Ltd (ScienceLab)**, Ambala, among its chemistry lab equipment in India.

Specifications

Parameter	Detail
Model No.	SCL-581/040
Material	Ceramic
Surface	Glazed; crisp and smooth
Colour	White
Shape	Square
Care	Easy to clean and maintain
Tile size	Not stated; confirmed on enquiry

Typical Use

In acid–base titrations with indicators such as methyl orange or phenolphthalein, place the tile on the base of the burette stand directly under the flask. Against the white glaze the first permanent tint is much easier to see, so students are less likely to overshoot the end point. In redox titrations that use an external indicator – for example ferrous ammonium sulphate against potassium dichromate with potassium ferricyanide – drops of indicator are placed in a row on the tile, and a drop of the titration mixture is touched to one of them after each addition until the colour no longer appears.

In a volumetric analysis kit for a class, allow one tile per burette so every student works against the same background. Rinse the tile with distilled water and wipe it dry after each titration so that no trace of indicator is left to confuse the next reading. The tile is not a heating surface: glazed ceramic can crack if it is heated suddenly, so keep it off burners and hot plates and use a wire gauze on a tripod for heating.

Supplier

ScienceLab is the trading name of Jain Scientific Equipments Private Limited in Ambala, Haryana. The company supplies accessories like these glazed tiles, so a titration station can be equipped completely by institutions in India and international markets. The [company profile](#) has more information.

Rela

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

Thes

1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India



Email : sales@scilabexport.com

Phone: +91-7082934803

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India