



Email : sales@scilabexport.com

Phone: +91-7082934803


Product Name :

Ferrous Ammonium sulphate

Product Code :

SCL-847/167

Image not found or type unknown

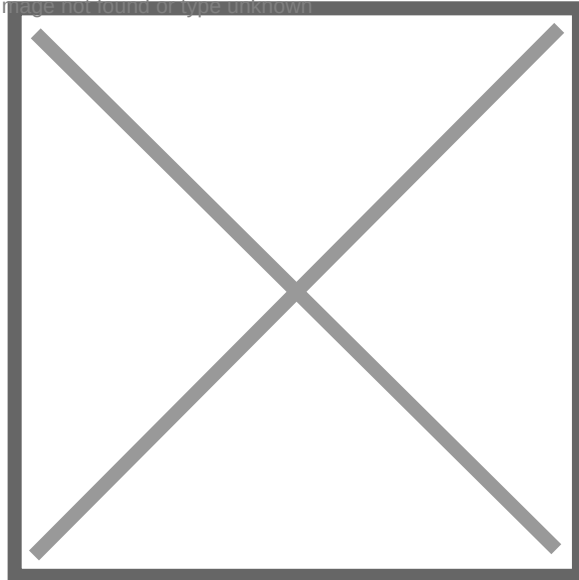

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India



Email : sales@scilabexport.com

Phone: +91-7082934803

Description :

Ferrous Ammonium sulphate (Product Code: SCL-847/167), widely known as Mohr's salt, is a pale green iron(II) double salt used to prepare standard solutions for permanganate and dichromate titrations in chemistry laboratories. It resists air oxidation better than plain ferrous sulphate, which is why senior secondary schools and colleges weigh it out directly; each pack holds 500 gm of analytical reagent with an MSDS. It is supplied by **Jain Scientific Equipments Pvt Ltd (ScienceLab)**, Ambala, India, a laboratory chemicals supplier in India.

Specifications

Parameter	Detail
Product code	SCL-847/167
Common name	Mohr's salt (ferrous ammonium sulphate)
Grade	Analytical reagent
Pack	500 gm
Safety data sheet	Supplied (MSDS)
Shelf life	Shelf life depends on the chemical, its grade and storage conditions; ask us for the expiry date of the batch you are quoted.
Supplied by	Jain Scientific Equipments Pvt Ltd (ScienceLab), Ambala, India

What Is Ferrous Ammonium Sulphate (Mohr's Salt) Used For?

The best-known practical is the determination of the strength of a potassium permanganate solution. Students weigh Mohr's salt accurately on a chemical balance, dissolve it in water containing some dilute sulphuric acid, make it up to the mark in a volumetric flask, and titrate portions against permanganate from a burette. No separate indicator is needed: the first lasting pale pink colour marks the end point. This experiment, together with its oxalic acid counterpart, forms the core of the redox section of a typical titration kit for schools.

Dilute sulphuric acid is added when dissolving the salt to stop hydrolysis and slow oxidation of iron(II) by air, so make the solution on the day of use. The salt also features in crystallisation practicals, where students prepare Mohr's salt from ferrous sulphate and ammonium sulphate. Keep the stock container closed and dry; crystals that have turned yellowish-brown on the surface have partly oxidised and are no longer suitable as a standard.

Supplier

Jain Scientific Equipments Private Limited, which presents itself to customers as ScienceLab, supplies Ferrous Ammonium sulphate from Ambala, Haryana. The company sources titration reagents such as Mohr's salt and quotes them with burettes, flasks and other glassware for institutions in India and international markets. See the [company profile](#).

Related Chemistry Lab Products

Reagent **Website:** <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

• 1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt, Harviana 133001, India. [Contact Us](#)



Email : sales@scilabexport.com

Phone: +91-7082934803

Website: <https://www.scilabexport.com>, **Email:** sales@scilabexport.com

1226/1-5, Bengali Mohalla, Science Market, Ambala Cantt. Harvana 133001, India